


SILEX[®]

Sacroiliac Joint
Fusion System



- ❑ 12.5mm Anchor Implants offered (single-thread and dual-thread)
- ❑ 7mm Locking Implants offered (dual-thread)
- ❑ 12.5mm Anchor Implants are fenestrated to allow for direct bone graft apposition
- ❑ Decortication and direct visualization of the Sacroiliac Joint achieved through innovative Silex[®] procedure

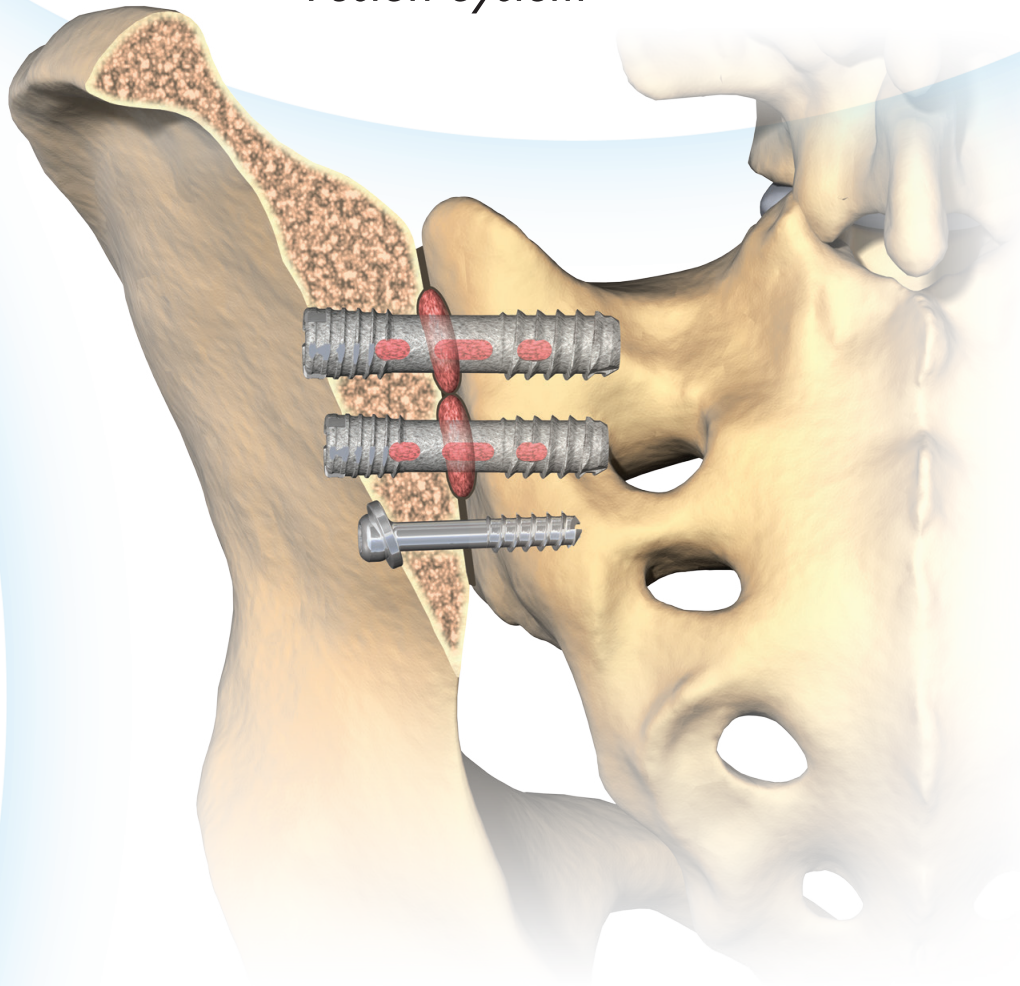



12.5mm Anchor Implants

Available in 30mm - 70mm lengths (5mm increments)

- ❑ Single-thread and dual-thread implants
- ❑ All implants titanium plasma-coated to aid in frictional resistance and provide bioactive template
- ❑ All implants fenestrated to allow for direct bone graft apposition
- ❑ Multiple lengths allow for anatomical variation
- ❑ Dual-thread to aid in compression across the joint
- ❑ Anti-migratory scallop features to help prevent implant back out
- ❑ All implants are self tapping
- ❑ All implants cannulated



7mm Locking Implants

Available in 30mm - 70mm lengths (5mm increments)

- ❑ Dual-thread implant only
- ❑ Optional 13mm washer available
- ❑ Multiple lengths allow for anatomical variation
- ❑ Anti-migratory scallop features to help prevent implant back out
- ❑ All implants are self tapping
- ❑ All implants cannulated



13mm Washer

- ❑ Allows for distribution of load
- ❑ For use with 7mm Locking implant



X-spine Systems, Inc.


452 Alexandersville Rd., Miamisburg, OH 45342
Phone: 800-903-0640 • Direct: 937-847-8400 • Fax: 937-847-8410
www.x-spine.com

INDICATIONS: See Package Insert for a more complete listing of indications, contraindications, warnings, precautions, and other important information.

LIMITED WARRANTY and DISCLAIMER: X-spine products have a limited warranty against defects and workmanship and materials. Any other express or implied warranties, including warranties of merchantability or fitness, are disclaimed.

WARNING: In the USA, this product has labeling limitations. See package insert for complete information.

CAUTION: USA Law restricts these devices to sale by or on the order of a physician.

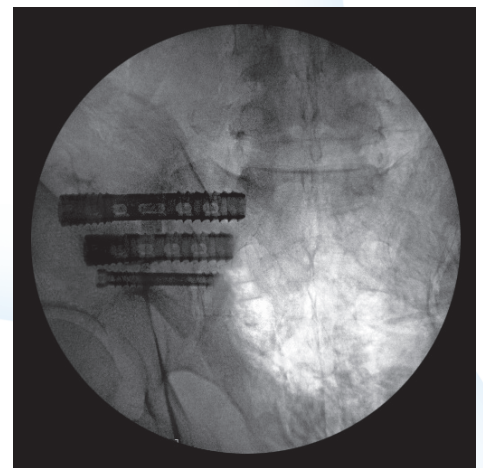

X-spineSM the X-spine logo and Silex[®] are trademarks or servicemarks of X-spine Systems, Inc.

Patents Pending.

All products are not currently available in all markets.
© 2015 X-spine Systems, Inc., All rights reserved.

CE0086

F-1000.24 Rev C 04/15

